



# SeekMate Tinitan<sup>TM</sup>

## Fluorescence Cell Counter



**A Smart Tiny Titan Is  
Counting Cells**


---

Compact size  
Intelligent detection



## SeekMate Tinitan™

SeekMate Tinitan™ FL is a highly automated fluorescence cell counter developed independently by SeekGene Biosciences. It integrates both bright field and fluorescence modes, offering a cutting-edge solution for high-efficiency fluorescence cell counting. Equipped with the latest optical components and a dual focusing mode, it captures super clear images. The incorporation of AI deep learning and personalized parameter adjustment provide accurate and stable cell counting solutions for various cell types (including primary cells and cell nuclei). The device is compact and compatible with mobile power sources, making it convenient to carry outdoors for cell counting work in any laboratory setting. SeekMate Tinitan™ FL is a handheld "counting expert" with precise counting, full function, flexibility and convenience.

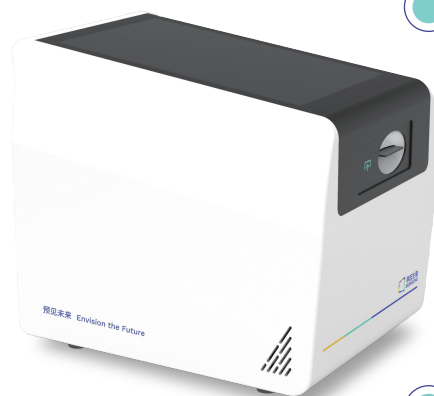

Compact size

Dual focusing mode

Well-trained AI

Sample-saving

Type-C interface

## Dual focusing mode

The dual focusing mode combines automatic and manual focusing, specifically designed for primary cells and special samples. By presetting optimal focusing parameters, it achieves automated accurate counting, overcoming challenges like diverse cell types, varying sizes, irregular shapes, poor compatibility with fluorescent dyes, and interference from impurities and debris.



Brightfield and Fluorescence image view of K562 cells


## More diverse cell types

With its broad applicability to cell sizes ranging from **0.5 to 180  $\mu\text{m}$** , SeekMate Tinitan™ FL has been successfully validated across diverse cell cultures, including primary cell lines and PBMCs.



## Well-trained AI algorithm

Utilizing a database of nearly **10,000** cell images to train dozens of cell detection parameters with one-touch access, it provides rapid detection and intelligent classification of cell aggregation, impurities, fragments, and cells with no nucleus.

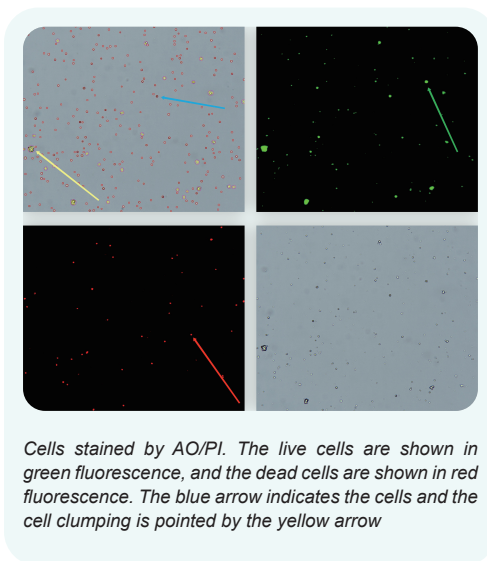


## Compact all-in-one design

The all-in-one machine package eliminates the need for an external computer, saving experimental space. With a compact body of **22.6\*15.8\*19.9 cm**, it is more convenient to carry outdoors. The 8-inch touchscreen provides a beautiful interface with clear graphics and flexible operation.

## User-friendly Type-C power interface

The Type-C power interface is compatible with power bank (25V 3.25A), extending continuous outdoor operation.

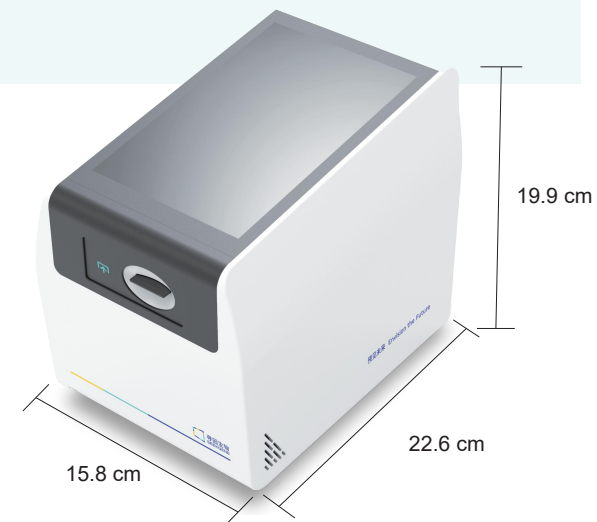


The compact size of Tinitan saves your lab bench space






## Accurate and stable counting


Validation studies using the SeekMate Tinitan™ FL and a hemocytometer have shown excellent accuracy and repeatability across various cell types, with CV values  $\leq 5\%$  and a linear correlation coefficient exceeding 99%.

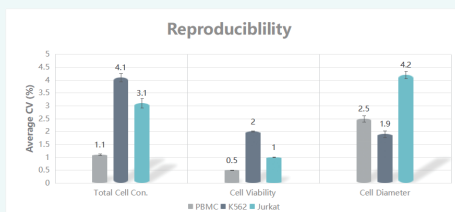

Great reproducibility was indicated for Tinitan through 10 repeated testings using PBMC, K562 and Jurkat cells, with the CV value of lower than 5%

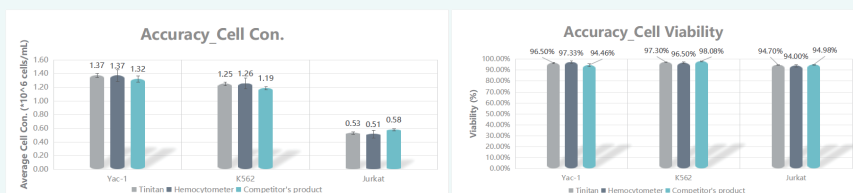

Testing of Tinitan, hemocytometer, and the competitor's product shows high comparability in cell concentration and viability.

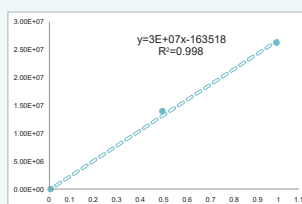


Dilution series of a K562 cell sample. The linear regression value suggests the high accuracy of Tinitan

## Multiple savings design

Small sample volume of **10  $\mu\text{L}$**  saves samples and reagents, reducing consumable costs. Tinitan can handle up to **6** samples simultaneously, automating processes and minimizing manual labor.



10  $\mu\text{L}$  sample loading to the chip

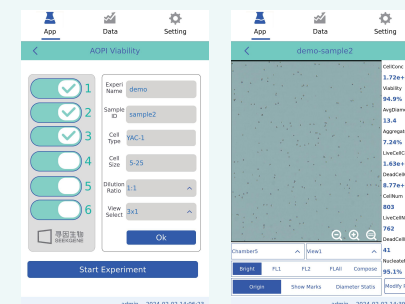


Setting parameters before detection

Result and image overview of YAC-1 cell line

## Compliance with FDA 21 CFR Part11 regulations

Compliance with FDA 21 CFR Part11 regulations ensures adherence to stringent standards, encompassing multi-level permission management, event tracking, signatures, report generation, and electronic record archiving, guaranteeing data security and integrity.

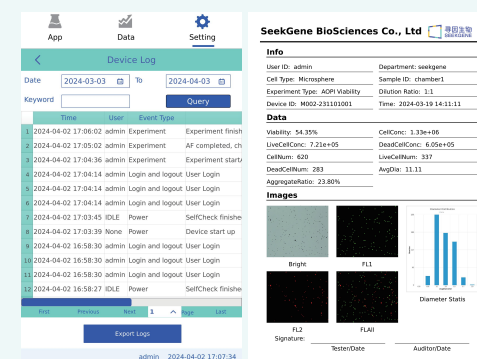

# Technical Specifications

| Parameters                       | Description                                                                                                                                          |
|----------------------------------|------------------------------------------------------------------------------------------------------------------------------------------------------|
| Dimensions                       | 22.6 x 15.8 x 19.9 cm                                                                                                                                |
| Weight                           | 4.5 kg                                                                                                                                               |
| Total cell concentration range   | 1x10 <sup>4</sup> - 3x10 <sup>7</sup> cells/mL                                                                                                       |
| Optimum cell concentration range | 5x10 <sup>4</sup> - 1x10 <sup>7</sup> cells/mL                                                                                                       |
| Cell size range                  | 0.5 - 180 µm                                                                                                                                         |
| Cell types                       | Cell lines, primary cells, PBMC, small cells, cell nuclei, spheroids, and other cell types                                                           |
| Count mode                       | Bright-field and Fluorescent                                                                                                                         |
| Sample volume                    | 10 µL                                                                                                                                                |
| Counting time                    | Trypan blue 25 s; AO/PI 40 s                                                                                                                         |
| Counting fields                  | 3                                                                                                                                                    |
| Sample throughput                | 6                                                                                                                                                    |
| Fluorescence channels            | Ex/Em: 470 nm/535 nm<br>Ex/Em:470 nm/600 LP<br>Ex/Em: 525 nm/600 LP                                                                                  |
| Camera                           | 5 mega pixels                                                                                                                                        |
| Focus method                     | Auto + Manual                                                                                                                                        |
| Detected indicator               | Cell count, cell viability, cell concentration, nucleation rate, clumping rate, mean cell diameter, transfection rate, growth curve, data comparison |
| Touch screen                     | 8 inch HD resolution                                                                                                                                 |
| Data storage                     | 500 G                                                                                                                                                |
| Data export method               | 1 x USB 3.0, Wi-Fi, cloud server                                                                                                                     |
| Power supply                     | Type-C power interface; mobile power with built-in battery                                                                                           |
| Objective magnification          | 5X                                                                                                                                                   |

# Ordering Information

| Name                                        | Quantity     | Art. No. |
|---------------------------------------------|--------------|----------|
| SeekMate Tinitan™ Fluorescence Cell Counter | 1 Instrument | M002C    |
| Cell Counter Chip                           | 50 chips/Kit | SP00022  |
| AO/PI Staining Kit                          | 5 mL/Kit     | K01701   |



Follow our social media to stay ahead with latest advances in single-cell field



Mail: [info@seekgene.com](mailto:info@seekgene.com)

Website: <https://www.seekgene.com>